

**AQUATEC®**

# The MAJOR

**Aquatec** has been providing safe bathing options for over 2 generations of users. Never stoop to lift someone out of the tub again. The **Aquatec Major** will rise to meet you at the top of the tub, then slowly and gently lower you into the warm bath waters. The **Aquatec Major** is a full-length bathlift that provides a bathing platform for youth or adult who cannot be placed in a seated position. It is also ideal for group facilities where bathers have a variety of disabilities and positioning needs. With a top capacity of 500 pounds, the **Major** can be used for bariatric clients.



## *Components*

The **Major** is constructed of stainless steel and plastic components. The bathlift is custom built to the dimensions of your tub. The unit includes all parts necessary for basic operation, base, top platform with cover mat, hand control, water supply hose and quick disconnect fitting.

The lift consists of a lifting base with 2 lifting scissors and bladders. The top panel with cover mat has 4 folding edges that will open at the top of your tub to form a solid platform for towel drying or dressing. The blue cover mat is non-latex mesh, made from blown polyvinyl chloride (PVC). It is removable and machine washable. These units are custom built and take several weeks for delivery. Your sales representative can assist with the details and order forms needed for ordering a **Major**.

## *Operation*

The **Aquatec Major** is operated with a joystick style hand control. The optional reclining back supports are also controlled from the hand control. The **Major** will lower to 2.75" from the bottom of a standard tub.

## *Warranty*

**Aquatec** offers a limited 2 year warranty on the **Major** to be free of manufacturer's defects. This warranty excludes covermats.

**CLARKE HEALTH CARE PRODUCTS, INC.**

1003 International Drive • Oakdale, PA 15071

Phone (888) 347-4537 • Fax (724) 695-2922 • [www.clarkehealthcare.com](http://www.clarkehealthcare.com)

The MAJOR

**AQUATEC®**

# The MAJOR

## Options/Accessories

The **Major** can be ordered with any of 4 back support options: stationary or reclining, with or without adjustable side laterals. A seat belt can be mounted, but due to safety concerns you cannot order a chest belt without a back support. Height adapters can be used for tubs deeper than 19" for a maximum height of 22". All the above options must be ordered with the original lift.



Straight Back



Reclining Straight Back



Head Support and Belts



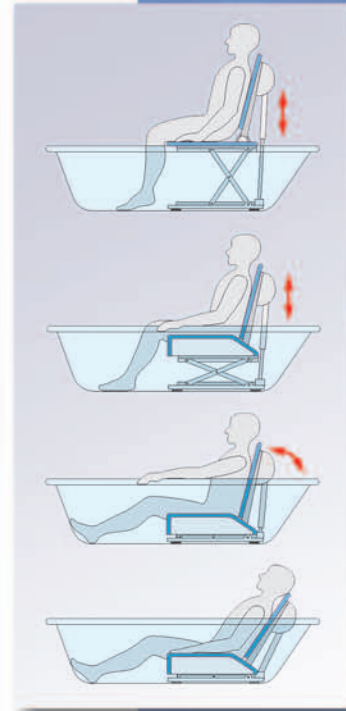
Special Back



Reclining Special Back



Sand Wedge Cushion

**SPECIFICATIONS****MAJOR**

Length of Seat	Custom
Width of Seat	Custom
Footprint, at bottom of tub	Custom
Width of top platform, opened	Custom
Fits tub with maximum width at top	40"
Maximum height of lift	19"
Height with Adapters (1,2, or 3")	20-22"
Fits tub with maximum depth (lift with height adapters)	22"
Height of back support, all 4 styles	16"
Width of straight back support	12.5"
Width of back support with laterals, flaps open	28.5"
Lifting Capacity	550 lbs.
Weight of Lift	NA
Product Number	A152112



Special Cushion



Rotary Seat with Transfer Board

**CLARKE HEALTH CARE PRODUCTS, INC.**

1003 International Drive • Oakdale, PA 15071

Phone (888) 347-4537 • Fax (724) 695-2922 • [www.clarkehealthcare.com](http://www.clarkehealthcare.com)

# The MAJOR