



The SYMPHONY

Symphony is a elegant orchestration of quality and durability and is the finest rollator on the market today. The handles have stepless adjustment and anatomic hand grips. These grips require less grip strength and reduce hand fatigue especially for arthritic hands. **Symphony** has a flip-up seat for greater reach or for gait training.

Frame

The lightweight aluminum frame is welded for greater durability. No screws to loosen at the joints. **Symphony** is the strongest aluminum frame available today. Grey anodized finish. **Symphony** is available in 3 frame heights. Easy to fold frame has a safety latch to prevent accidental folding. Aluminum frames are suitable for corrosive environments, such as coastal areas. The frame is warranted for 2 years.

Braking

Symphony comes with the patented Dolomite Automatic Braking System with adjustable tension for varied grip strength. Brake handles have wide flange and large loops for comfortable hand positioning. When parked and using the seat, push down the handles to set the parking brakes. The brake mechanism is located within the fork for protection. Accessory slow down brake provides adjustable tension for hilly terrain or to even out the direction of travel for stroke patients. The brake handle assembly is warranted for 5 years.

Standard Features


- Flip-up seat for gait training and greater reach
- Curb climber on rear wheel assists when lifting front of walker
- Anatomic grips reduce hand fatigue; can be rotated inward for a more natural position
- Adjustable back support
- Basket for shopping
- 8" indoor/outdoor wheels
- Stepless handle height adjustment, adjustment marks printed on handle
- Hidden brake cables



CLARKE HEALTH CARE PRODUCTS, INC.

1003 International Drive • Oakdale, PA 15071

Phone (888) 347-4537 • Fax (724) 695-2922 • www.clarkehealthcare.com

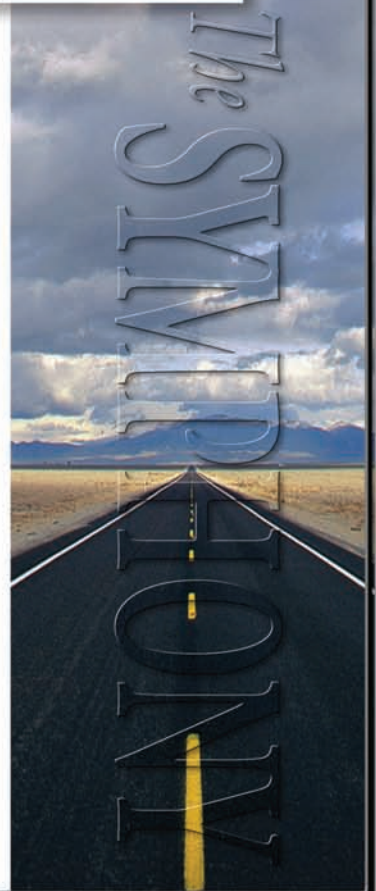



The SYMPHONY

Accessories



OXYGEN TANK HOLDER



ONE-HAND BRAKE



SLOW DOWN BRAKE



FOOD & BEVERAGE TRAY



CANE HOLDER



HEMIPLEGIC HANDLE



UNDER SEAT BAG

Accessories can enhance the use of your walker – the Tray gives you a platform to carry food, the Oxygen Tank Holder will safely support a B, C or D tank. A Cane Holder can secure a cane or crutch to the side of the walker, Slow Down Brakes have a tension control to allow individual settings for each wheel. One-Hand Brakes can control your braking from one handle. The Hemiplegic Handle allows you to use the body to push for momentum while braking is done with one hand. Under Seat Bag gives you a place to stash items out of sight.

SPECIFICATIONS	SUPER LOW	LOW	STANDARD
Dist. Between Rear Wheels	18.5"	18.5"	18.5"
Minimum Handle Height	25"	27.5"	30"
Maximum Handle Height	30"	33"	36"
Seat Height	18"	21"	24"
Overall Width	23"	23"	23"
Overall Length	22"	24"	27"
Width Between Handles	17"	17"	17"
Walker Weight (lbs.)	15	15	15
Product Number	D12162	D12161	D12160

CLARKE HEALTH CARE PRODUCTS, INC.

1003 International Drive • Oakdale, PA 15071

Phone (888) 347-4537 • Fax (724) 695-2922 • www.clarkehealthcare.com

