

AQUATEC® *The AQUATEC Bathlift*

*Aquatec shown with
Rotary Seat with
Transfer Board*



The latest technology for safe and carefree bathing. The **Aquatec** bathlift comes from a family of high quality bath products.

Aquatec is a battery powered bath seat that will lower and raise you again in your tub. Enjoy a long soak and the therapeutic effects of a warm bath. With a reclining back, optional seat widths and back laterals, the **Aquatec** is perfect for those with special bathing needs. The hand remote has easily identified up and down buttons and low power indicator. The lithium-ion battery is in the hand control for convenient recharging. Self release suction cups adhere securely to smooth tub floors and the seat flaps will rest on the top of the tub to allow a stable transfer.

Features

- Backrest with built-in neck support
- Lithium-ion battery with discharge protection, has no memory effect
- Self release suction cups
- Reclining back support up to 40°

Non-slip surface



Components

Aquatec models include a unit to fit standard size tubs, **Aquatec Wide** model for whirlpool style tubs, a back with adjustable side laterals, known as a RSB model and **AquatecXL** model with a 400 pound weight capacity.

All models use a state-of-the-art hand control with integrated battery. The latest lithium-ion manganese-based rechargeable battery pack has total discharge protection and no memory effect. A low battery indicator will light on the hand control when the battery needs to be recharged. A safety feature will not allow the lift to lower when there is not enough power to rise again.

A high backrest with molded neck support can be folded forward by just one hand. The maintenance-free actuator is both powerful and silent. **Aquatec** has a stable and rust-free frame, reinforced lifting scissors and bottom plate.

Aquatec has a non-slip textured surface and comes with removable, washable covers. The bathlift disassembles into easy to carry parts; back, seat and hand control.

AquatecRSB




CLARKE HEALTH CARE PRODUCTS, INC.

7830 Steubenville Pike • Oakdale, PA 15071

Phone (888) 347-4537 • Fax (724) 695-2922 • www.clarkehealthcare.com

The AQUATEC Bathlift Accessories

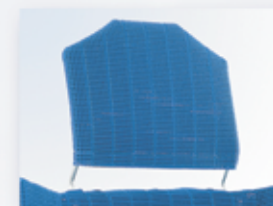


A15461 / A15462
Head support, height adjustable mounts on any back support



A1471210
Seat belt, adjustable cloth strap attaches to seat base


A1471209
Chest belt, used with side laterals, adjustable strap attaches to back support

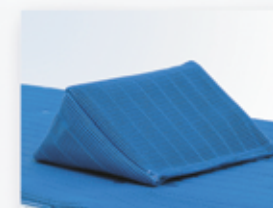

A13693 / A13694
Sand wedge cushion will stop forward slide associated with high tone or spastic bathers



A14895 / A14901
Abductor block, mounts to lift base

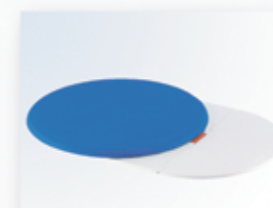

A403002
Rotary seat with transfer board, assists in positioning to center of lift seat. Center disc pivots and slides on base plate



A14586
Carry bag, water resistant nylon with carry handles, holds lift base and back support

AQUATEC®

The AQUATEC Bathlift Accessories



A16240
5" Suction pad provides greater area to grab on slightly textured surfaces



A15170 / A14694
Height adapters, adds up to 3 additional inches to base plate, order in 1" increments


Replacement parts including cover mats are also available.

Specifications:


DESCRIPTION	DIMENSIONS	WT. LBS.	THICKNESS
Head Support, Battery Powered Units	12"W x 8"H	1.5	1"
Sand wedge cushion	14"L x 9.5"W x 6"H	16	
Rotary seat with transfer	21"L x 12.5"W	8	.75"
Carry Bag	29"L x 12"W x 14"H	3	
Abductor, A14895	3"Dia. x 4.5"H	.6	
Abductor, A14901	5"Dia. x 5"H	1.6	

AQUATEC

AQUATEC[®] The AQUATEC Bathlift

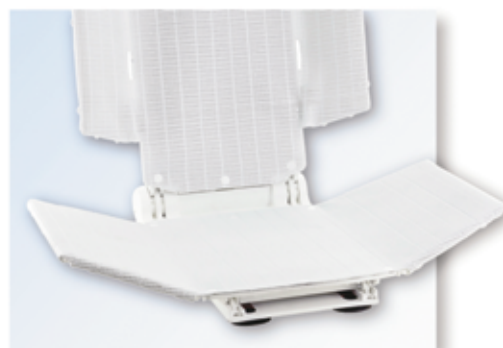

AquatecRSB Wide

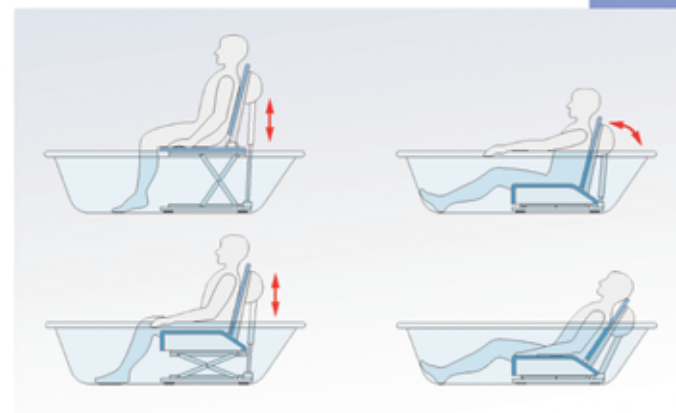

Warranty

Aquatec has a limited three year warranty against manufacturer's defects. Covers, battery and suction cups are excluded.

Specifications:

DESCRIPTION	SPECIFICATION				
	Aquatec	Aquatec Wide	Aquatec RSB	Aquatec RSB Wide	Aquatec XL
Seat width without sideflaps	14.75"	14.75"	14.75"	14.75"	14.75"
Seat width with sideflaps	27.88"	35.5"	27.88"	35.5"	27.88"
Sitting area length	19"	19"	19"	19"	19"
Maximum seat height	17"	17"	17"	17"	17"
Minimum seat height	2.5"	2.5"	2.5"	2.5"	2.5"
Backrest width	14"	14"	27"	27"	14"
Backrest height	26"	26"	26"	26"	26"
Total length (backrest upright)	24"	24"	24"	24"	24"
Total length (backrest reclined)	33"	33"	33"	33"	33"
Total height including backrest	42"	42"	42"	42"	42"
Bottom plate, footprint	12" x 23"	12" x 23"	12" x 23"	12" x 23"	12" x 23"
Cover thickness	.11"	.11"	.11"	.11"	.11"
Total weight without hand control	27 lb	29 lb	31.25 lb	31.25 lb	27 lb
Seat weight	17 lb	19 lb	17 lb	19 lb	17 lb
Backrest weight	10 lb	10 lb	12.5 lb	12.5 lb	10 lb
Hand control weight	1 lb	1 lb	1 lb	1 lb	1 lb
Load capacity	300 lb	300 lb	300 lb	300 lb	400 lb
Backrest reclines to	40°	40°	40°	40°	40°
Part Numbers Blue	A1471565	A1471651	A1471575	A1471582	N/A
Part Numbers White	A1471197	A1471650	A1471202	A1471579	A1471199

CLARKE HEALTH CARE PRODUCTS, INC.

7830 Steubenville Pike • Oakdale, PA 15071
Phone (888) 347-4537 • Fax (724) 695-2922 • www.clarkehealthcare.com